VELFORM VACU CLEAR™


USER MANUAL

EN.....	2	LV.....	26
AL.....	4	MD.....	27
BG.....	6	ME.....	29
BIH.....	10	MK.....	31
CZ.....	12	PL.....	33
EE.....	14	RO.....	35
HR.....	16	RS.....	37
HU.....	18	RU.....	39
KS.....	20	SI.....	41
KZ.....	22	SK.....	43
LT.....	24	UA.....	45

(EN) IMAGES / (AL) FOTOGRAFIN / (BG) ИЛЮСТРАЦИИ / (BIH) SLIKE / (CZ) OBRÁZKY / (EE) JOONISED / (HR) SLIKE / (HU) KÉPEK / (KS) FOTOGRAFIN / (KZ) СУРЕТТЕР / (LT) PAVEIKSLIAUKAI / (LV) ATTĒLI / (MD) IMAGINI / (ME) SLIKE / (MK) СЛИКИ / (PL) RYSUNKI / (RO) IMAGINI / (RS) SLIKE / (RU) RYSUNKI / (SI) SLIKE / (SK) OBRÁZKY / (UA) МАЛЮНКИ

Picture 1

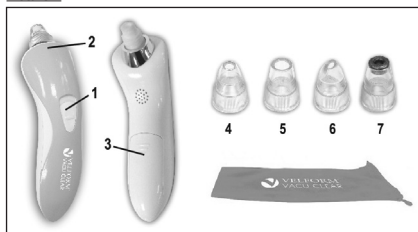
Picture 2

Picture 3

Picture 4

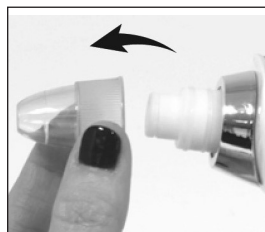Picture 5

Picture 6

Picture 7

EN

Please read this manual thoroughly before using and save it for future reference.

WARNINGS

1. Do not wear necklaces, earrings or any facial jewelry while using this product. Necklaces, earrings and other types of jewelry should be removed before using the product.
2. This product is not for use on areas of the body that have undergone any sort of surgery or have any scars.
3. If you experience pain while using this product, discontinue use immediately.
4. This is not a medical device. Nothing contained should be considered medical advice, diagnosis or treatment.
5. This product is not waterproof and should be not used in the shower or any situation that involves water.
6. Do not store this product in the bathroom or a room that has high moisture content.

7. Do not store or use this product near sources of heat or fire.
8. Keep out of reach of children. This device is not a toy and should not be used by children or stored accessible to them.
9. Do not use on irritated skin. Discontinue use if skin irritation occurs.
10. Exercise caution when using around eyes.
11. This device is intended for household and dry use only.
12. Do not submerge unit under water.
13. Never use the device without one of the suction heads correctly installed on the suction nozzle.
14. Regularly remove blockages caused by dust, fluff, hair etc. because they greatly reduce device's suction power and/or general performance.
15. Store in a cool dry place when not in use.
16. If the device is not working as it should, has received a sharp blow, has been dropped, damaged, left outdoors or dropped into water, do not use it.
17. This device should never be switched on when not in use, while replacing accessories or during cleaning.
18. In case of malfunction and/or damage do not attempt any repair by yourself and ensure that the product is properly disposed in accordance with local laws.
19. Although the device has been checked its usage and consequences are strictly user responsibility.

ABOUT THE PRODUCT

Wellneo® Velform Vacu Clear™ is fast and easy way to unclog and clean your pores without squeezing. The powerful vacuum sucks out the excessive matter resulting in removing whiteheads and blackheads. Wellneo® Velform Vacu Clear™'s powerful yet gentle vacuum technology penetrates deep leaving your skin clean and beautiful. Months and months of walking around with blackheads on your face or body are gone! Don't walk around with them - leave them to Wellneo® Velform Vacu Clear™.

SET INCLUDES (See Picture 1)

- 1x Main unit** (1 = power and setting switch; 2 = Suction nozzle; 3 = Battery compartment)
- 1x Small head** (4; gentle suction, perfect for everyday use. Good for use on sensitive skin and nose area.)
- 1x Large head** (5; strong suction, perfect for blackheads and deep pore cleaning. NOTE: Do not use around the eyes.)
- 1x Oval head** (6; oval opening for getting into corners around the nose and mouth, and for sensitive skin around the eyes.)
- 1x Peeling head** (7; Helps exfoliate dead skin while pulling dirt and oil from the skin.)
- 1x Travel microfiber pouch**

BATTERY INSTALLATION

- Step 1:** Turn the device off by sliding the power and setting switch all the way down.
- Step 2:** Open the battery compartment cover by sliding it in the direction of the arrow (Picture 2).
- Step 3:** Insert 2 AA batteries (not included) into the battery compartment. Make sure to line up the (+) and (-) markings on the battery with the (+) and (-) markings on the inside of the battery compartment (Picture 3).
- Step 4:** Reinstall (slide back) the battery compartment cover.

CHANGING SUCTION HEADS

- Step 1:** Turn the device off by sliding the power and setting switch all the way down.
- Step 2:** Place one hand on the suction head and the other hand on the main unit before firmly pulling the suction head from the suction nozzle (Picture 4).
- Step 3:** Firmly pull off the suction head and reinstall your desired suction head in the same manner (Picture 5).

USING

- Step 1:** Start by washing your face with warm water, steam or a hot wet towel to open up your pores.
 - Step 2:** Make sure to dry your face and hands completely before moving on to the next step.
 - Step 3:** Select your desired suction head (suction head selection tips/suggestions are available in the «set includes» section). Wellneo® Velform Vacu Clear™ can be turned on by sliding power and setting switch up. The middle position is low suction and the top position is high suction (Picture 6). You can always shut off the device at any time by sliding the power and setting switch all the way down.
 - Step 4:** Apply the Wellneo® Velform Vacu Clear™ so the head comes gently in contact with your skin (Picture 7). Slowly move the device around your face, arms or chest, without leaving it in one spot for more than two seconds. For best results, focus on your problem areas.
- NOTE: DO NOT move the device back and forth on a specific area or use it on the same spot for more than 2 seconds.*

MAINTENANCE

After each use, clean the suction head by rinsing it with warm water. Allow the suction head to air dry. Do not attempt to clean the main unit by rinsing it with water because it is not water-resistant. To clean the main unit, wipe it down with a soft, damp cloth.

HELPFUL TIPS

- After using Wellneo® Velform Vacu Clear™, make sure to apply moisturizer to the area the product was used.
- Before and after use, do not use facial scrubs or other exfoliating products.
- Do not press down heavily during use of Wellneo® Velform Vacu Clear™. Its suction capabilities are enough to clean your pores efficiently.
- Slight redness is normal and will disappear after short time.

NOTE: Wellneo® Velform Vacu Clear™ can be used once a week on dry type skin. On oily and mixed type skin, it can be used twice a week. Do not use it for longer than five minutes a day.

TROUBLESHOOTING

- If the product does not work, check if the batteries are inserted correctly into the device. If it still does not work, check if the suction head is attached correctly or any of the parts needs cleaning.
- If unit begins to slow down, clean the suction head of residue and/or replace the batteries if needed.

WELLNEO VELFORM VACU CLEAR

ID: 110033461

Model no.: VCSVFCIND0202

Rated voltage: 3.0 V DC (2 x 1.5 V AA alkaline battery - not included) / Power: 3 W

ORT no. and Production date can be found on the product.

Wellneo trademark owner, Importer to EU, Distributor for Studio Moderna: Top Shop International SA, Via al Mulino 22, CH-6814 Cadempino, Switzerland.

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