concept

Osobní váha
Osobná váha
Waga osobista
Személmérleg
Ќермеґа сvari
Personal scale

Personenwaage
Pèse-personne
Bilancia personale
Balanza personal
Cântar persoane

VO3000

CZ

SK

PL

HU

LV

EN

DE

FR

IT

ES

RO

PODĚKOVÁNÍ

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ACKNOWLEDGEMENTS

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

Technical specifications	
Battery	2x 1,5 V AAA
Weighing range	2 - 180 kg
Step	100 g
Unit	kg, lb, st
Memory	10 users
User mode	Adult, child (10-85 years)
Measured body parameters	Weigh, body fat, body water, muscle mass, bone mass, BMI, BMR
Measurement range	Body fat: 0-46,4 % Body water: 37.8-66 % Muscle mass: 20-56 % Bone mass: 2-20 %

IMPORTANT SAFETY PRECAUTIONS

- Do not use the appliance in any other way than as described in this manual.
- Remove all the covering and marketing materials from the appliance before the first use.
- Observe the marked polarity of the batteries when inserting and replacing them.
- If you do not use the appliance for a long time, remove the battery from the appliance.
- Do not allow children or unauthorized people to handle the appliance, use it out of their reach.
- Persons with reduced mobility, reduced sensory perception, mental disadvantage or persons not familiar with the operation must use the appliance only under the supervision of responsible person, who is familiar with the instructions.
- Pay extra attention if the appliance is used near children.
- Do not use the appliance as a toy.
- Only use the appliance on a horizontal, stable surface.

- Before cleaning, installing and dismantling the appliance covers, switch off the appliance and remove the batteries!
- Do not use coarse and chemically aggressive substances to clean the appliance.
- If the segments or digits are displayed incorrectly, try to reset the appliance. Remove the batteries, wait a few seconds, and insert the battery again. If this does not work, contact your dealer or an authorized service centre.
- Do not use the appliance if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- The appliance is only suitable for home use and is not intended for commercial use.
- Keep the appliance away from heat sources such as radiators, ovens, etc. Protect it from direct sunlight and humidity.
- Do not place the appliance in a microwave or other oven.
- Do not touch the appliance with wet hands.
- The appliance is a fragile measuring device. Prevent shocks and falls.
- Do not repair the appliance yourself. Contact an authorized service centre.
- The data measured by this scale are for information purposes only. They are not a medical advice.
- Do not use this appliance if you have a pacemaker, ECG or any other internal medical device.
- Do not use if you are pregnant.
- Do not use on a slippery floor or if you have wet feet.
- The data accuracy may vary if you weigh yourself after exercise.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Cleaning and maintenance performed by the user must not be carried out by children unless they are at least 8 years old and under supervision. Children aged under 8 years of age must be kept away from the appliance and its power cord. Do not let children play with the appliance.



Damages from failure to follow the manufacturer's instructions are not covered by the warranty.

PRODUCT DESCRIPTION

- 1 Weighing area
- 2 Display

Control

- 3 Settings button
- 4 Move up
- 5 Move down

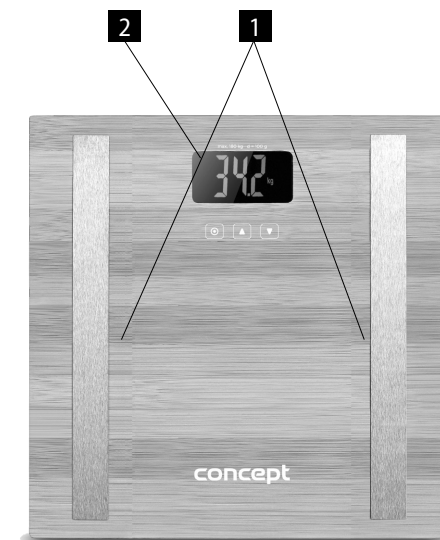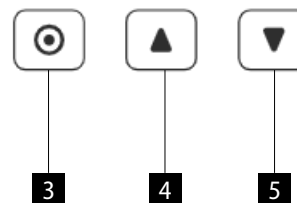OPERATING INSTRUCTIONS

BATTERY INSTALLATION

- Open the battery compartment cover on the bottom of the appliance.
- Insert 2 AAA batteries into the battery compartment and close the cover properly.

WEIGHT UNIT SELECTION

- Step on the scale to switch it on and use the up and down arrows to select the weight unit (kg, lb, st).

PERSONAL DATA SETTING

User number

Press the settings button and use the arrows to select the user number (options P0-P9), press the settings button to confirm and move to the next step.

Gender

The male/female icon flashes on the display, use the arrows to select the desired gender and press the settings button to confirm.

Age

Age flashes on the display, use the arrows to select your age (10-85 years), press the settings button to confirm.

Height (cm)

Height flashes on the display, use the arrows to and set your height (100-240 cm), press the settings button to confirm.

WEIGHING

- Place the scale on a hard, even floor.
- Stand with your bare feet in the middle of the scale, as shown in the Fig. 1.
- Your weight will appear on the display and flash.
- The scale switches off automatically after a few seconds.
- Use the "UNIT" switch on the bottom side of the scale to change the weight unit (kg, lb, st).

Fig. 1

BODY PARAMETERS MEASUREMENT

- Set personal data according to the instructions above.
- Wait until 0.0 kg (lb, st) appears on the display.
- Stand with your bare feet on the metal bars and wait, the body parameters are being measured.
- Body parameters gradually start to appear on the display.
- The parameters are measured using bioelectrical impedance analysis (BIA).

Comparison of BF (body fat) and BW (body water)

Age	Woman		Man		Note
	Fat %	Water %	Fat %	Water %	
≤30	4.0- 20.5	66.0- 54.7	4.0- 15.5	66.0- 58.1	Underweight
	20.6- 25.0	54.6- 51.6	15.6- 20.0	58.0- 55.0	Healthy weight
	25.1- 30.5	51.5- 47.8	20.1- 24.5	54.9- 51.9	Overweight
	30.6- 80.0	47.7- 13.8	24.6- 80.0	51.8- 13.80	Obesity
>30	4.0- 25.0	66.0- 51.6	4.0- 19.5	66.0- 55.4	Underweight
	25.1- 30.0	51.5- 48.1	19.1- 24.0	55.3- 52.3	Healthy weight
	30.1- 35.0	48.0- 44.7	24.1- 28.5	52.2- 49.2	Overweight
	35.1- 80.0	44.6- 13.8	28.6- 80.0	49.1- 13.80	Obesity

Comparison of the percentage of muscle and bone mass

Age	Woman		Man	
	Muscles (%)	Bones (%)	Muscles (%)	Bones (%)
16- 50	>34	>8	>40	>10
>50	>29	>8	>35	>10

Measured values are for information purposes only. If you have any questions concerning your health, consult your doctor.

DISPLAY

Err overload
 ----- underweight or overload
 LO replace the batteries

CLEANING AND MAINTENANCE

To clean the surface of the appliance, use only a damp cloth, no detergents or hard objects, as they can damage the surface of the appliance!

Never clean the appliance under running water, do not rinse it or immerse it in water!
Store the appliance in a dry place!

SERVICE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified specialist or the service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan:

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (Waste Electrical and Electronic Equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead, it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

Disposal of batteries and power supplies

If the appliance contains an accumulator or a battery cell, these components contain substances which may be harmful to the environment when disintegrated. Before disposing of the appliance or replacing the accumulator or the battery cell, it is necessary to remove the original component and hand over to a suitable collection point.

CAUTION!**Accumulator or battery cell:**

- Do not throw into a fire!**
- Never disassemble!**
- Never shorten the wires!**

If leakage occurs when the accumulator is damaged, the electrolyte is of an alkaline character. If the electrolyte gets on your skin, wash it gently with soap and water and seek medical attention.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without prior notice and we reserve the right to change them.

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Warranty

The manufacturer (or importer) is responsible for ensuring the product complies with the requirements of applicable legal regulations as well as those of the relevant technical standards. Moreover, they are responsible for ensuring the product has the properties the manufacturer described in documents related to the goods or those reasonably expected by the customer with regard to the nature of the goods or based on advertising produced by the manufacturer, and further they are responsible for ensuring the product is fit for the purpose proposed by the manufacturer or that a product of the same type is normally used for.

The quality warranty term is 24 months from product takeover by the customer.

The warranty does not apply to wear and tear caused by regular use. The customer shall not be entitled to any warranty claims if, prior to taking the product over, they knew the product contained a defect or if the defect is attributable to the customer.

The warranty specifically does not apply:

- if the product installation, operation and service conditions stipulated in the product operating manual have not been adhered to,
- to malfunctions caused due to mechanical, heat or chemical damage, short circuit, over voltage or incorrect installation,
- to malfunctions caused by an inexperienced third-party intervention,
- to malfunctions caused by natural disaster,
- to malfunctions caused by insufficient or inappropriate maintenance in violation of the operating manual, including malfunctions caused by water and other sediments,
- to changes in colour of the heating surface or to scratching of the surface caused as a result of using the products in an unusual manner,
- to appearance and functional changes caused by exposure to sunlight, thermal radiation of water and other sediments,
- if the service life of certain product parts expires, e.g. for accumulators, bulbs, etc.

The warranty does not apply to any products and services provided along with the product (gifts, promotional articles, etc.).

Filing a complaint

A complaint against a product defect must be filed as soon as identified, yet no later than before the end of the warranty term.

The customer must file a product complaint with the dealer from which they have purchased the product, or with any authorised service centre, a list of which is included in the product package, or available at www.my-concept.com.

While filing a product complaint, the product must be duly cleaned and securely packed so as to prevent any damage during its transport to an authorised service centre, where relevant, unless the product is delivered in person.

The customer must submit proof of having concluded a purchase contract for the product by producing the receipt.

While filing their complaint, the customer must indicate the noted defect and identify the preferred complaint application method.

Complaints processing

As long as the noted defect may be removed, the user has the right to have the defect duly removed free of charge on a timely basis.

Where such a procedure is not reasonable with regard to the nature of the defect, the user may require to be supplied a new defect-free product (replacement), or, where the defect applies to a part of the product only, replacement of the part concerned. However, if replacement of the product or any part thereof is not proportionate with regard to the nature of the defect, especially if the defect can be removed without undue delay, the customer has the right to have the defect removed free of charge.

If the noted defect is not removable, or if the customer becomes entitled to replacement of the product or a part thereof, yet the replacement is not possible, for example due to the product having been sold out, the customer has the right to return the product (withdrawal from the contract).

The customer shall have the right to a new product (replacement) or to replacement of a part of the product even if the defect can be removed, provided they cannot properly use the product due to repeated occurrence of the defect or due to a high number of such defects. In such a case, the customer also has the right to return the product (by withdrawing from the contract).

If the product is not returned (the customer does not withdraw from the contract), or if the customer does not apply the right to a new defect-free product (replacement), to replacement of a part thereof or to repair of the product, they may request a reasonable discount. The customer also has a right to a reasonable discount if a new defect-free product cannot be supplied to them, or if a product part cannot be replaced or the product repaired unless the situation is remedied within a reasonable time limit, or if remedying the situation would create major discomfort on the part of the customer.

The seller, authorised service centre or a staff member authorised by them must decide about each complaint immediately or within three business days in complicated cases. This term does not include a reasonable period of time, depending on the type of product concerned,

required for the defect to be assessed by an expert. A complaint, including defect removal, must be dealt with without any undue delay, yet no later than within 30 calendar days of the complaint filing date, unless the seller and the customer agree on a later deadline.

On returning the product (withdrawing from the contract) the customer must return any accessories and documents supplied along with the product.

The customer does not have the right to keep the defective parts and components of the product replaced as part of a repair of the product.

This shall be without prejudice to any other rights the customer may have in relation to the purchase of the product.

Remark: Complaints against products damaged in transport are governed by the carrier's complaints procedure.

Manufacturer

Jindřich Valenta - Concept
Vysokomýtská 1800, 565 01 Choceň
Czech Republic
tel.: +420 465 471 400
fax +420 465 473 304
Company ID No. 13216660
email: servis@my-concept.cz
www: www.my-concept.com

Product data

Model:	
Production number:	
Date of purchase:	Seal and signature of vendor:

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Verners VT Ltd.
S.C. DANTE INTERNATIOANL S.A.