concept

Osobní váha

Osobná váha

Waga osobista

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Çermeņa sviri

Personal scale


Personenwaage

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Bilancia personale

Balanza personal

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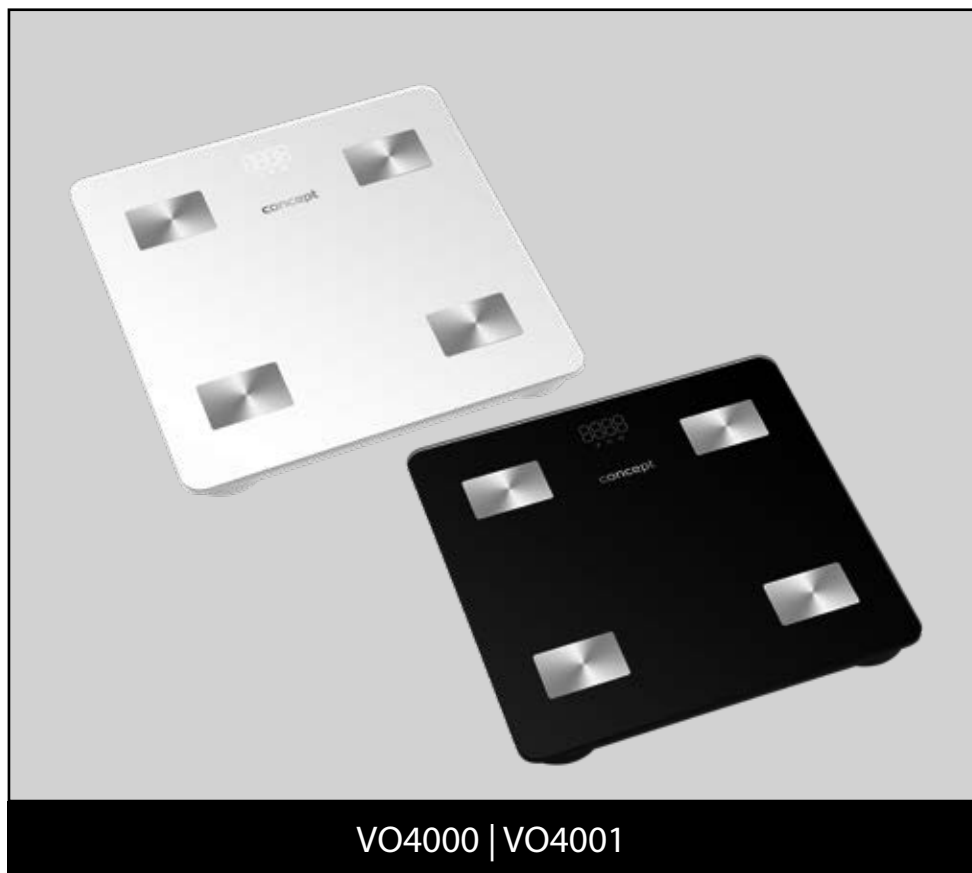VO4000 | VO4001

CZ

SK

PL

HU

LV

EN

DE

FR

IT

ES

RO

PODĚKOVÁNÍ

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ACKNOWLEDGEMENTS

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

Technical specifications	
Battery	3x 1,5 V AAA batteries (included in the package)
Weighing range	0.2 - 180 kg
Step	100 g
Frequency range	2.4 - 2.4835 GHz
Transmitter power	≤20 dBm

IMPORTANT SAFETY PRECAUTIONS

- Do not use the appliance in any other way than as described in this manual.
- Remove all the covering and marketing materials from the appliance before the first use.
- Observe the marked polarity of the batteries when inserting and replacing them.
- If you do not use the appliance for a long time, remove the battery from the appliance.
- Do not allow children or unauthorized people to handle the appliance, use it out of their reach.
- Persons with reduced mobility, reduced sensory perception, mental disadvantage or persons not familiar with the operation must use the appliance only under the supervision of responsible person, who is familiar with the instructions.
- Pay extra attention if the appliance is used near children.
- Do not use the appliance as a toy.
- Only use the appliance on a horizontal, stable surface.
- Before cleaning, installing and dismantling the appliance covers, switch off the appliance and remove the batteries!
- Do not use coarse and chemically aggressive substances to clean the appliance.
- If the segments or digits are displayed incorrectly, try to reset the appliance. Remove the batteries, wait a few seconds, and insert the


battery again. If this does not work, contact your dealer or an authorized service centre.

- Do not use the appliance if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- The appliance is only suitable for home use and is not intended for commercial use.
- Keep the appliance away from heat sources such as radiators, ovens, etc. Protect it from direct sunlight and humidity.
- Do not place the appliance in a microwave or other oven.
- Do not touch the appliance with wet hands.
- The appliance is a fragile measuring device. Prevent shocks and falls.
- Do not repair the appliance yourself. Contact an authorized service centre.
- The data measured by this scale are for information purposes only. They are not a medical advice.
- Do not use this appliance if you have a pacemaker, ECG or any other internal medical device.
- Do not use if you are pregnant.
- Do not use on a slippery floor or if you have wet feet.
- The data accuracy may vary if you weigh yourself after exercise.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Cleaning and maintenance performed by the user must not be carried out by children unless they are at least 8 years old and under supervision. Children aged under 8 years of age must be kept away from the appliance and its power cord. Do not let children play with the appliance.

Damages from failure to follow the manufacturer's instructions are not covered by the warranty.

PRODUCT DESCRIPTION

- 1 Weighing area
- 2 Display

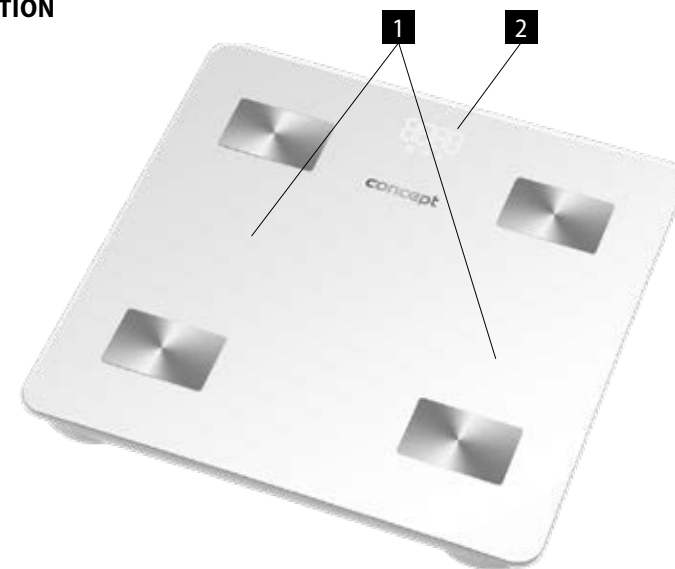OPERATING INSTRUCTIONS

BATTERY INSTALLATION

- Open the battery compartment cover on the bottom of the appliance.
- Insert 4 AAA batteries into the battery compartment and close the cover properly.

Note: After inserting the batteries, you can use the scale immediately. The scale displays only the current weight. Use the mobile app to measure other parameters.

DOWNLOADING THE APP

- Download the Feelfit mobile app to use all functions of the scale.
- You can download the app from Google Play (Android) or App Store (iOS). Scan one of the QR codes below to save time.

Google Play

App Store

PAIRING THE APP

- Enable the Bluetooth function in your mobile phone.
- Open the Feelfit app and follow the instructions on the screen.
- Select the basic data and register.
- Click „+“ on the display of the phone to add a new device and step on the scale. If the pairing is successful, an icon of a scale appears on the display of the phone, click on it.

**MEASUREMENT**


- Place the scale on a hard, even floor.
- Stand with your bare feet in the middle of the scale, as shown in the picture below.
- Wait for your weight to display. You will see a moving zero, which signals the measurement of body parameters. After the measurement is finished, your weight is displayed again.
- All measured data are displayed in the app.

**TIPS FOR OPTIMAL MEASUREMENT**

- To ensure the most accurate and consistent measurement, we recommend using the scale at the same time every day.
- Place the scale on a firm, dry floor. Weigh barefoot and position your feet so that your weight is evenly distributed on the scale. Wash and dry your feet.
- Do not weigh after exercising, showering, bathing or sauna, drinking alcohol, if you have fever or if you have just eaten.

DISPLAY

- The picture below shows the individual displays of the appliance.

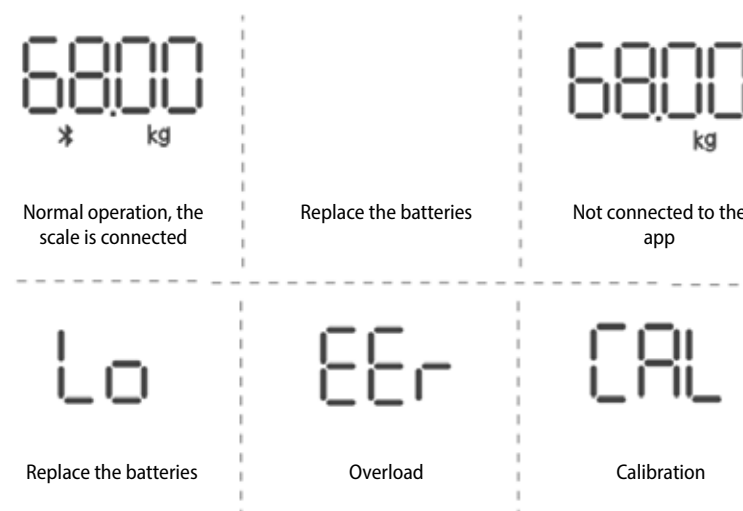**CLEANING AND MAINTENANCE**

To clean the surface of the appliance, use only a damp cloth, no detergents or hard objects, as they can damage the surface of the appliance!

Never clean the appliance under running water, do not rinse it or immerse it in water!
Store the appliance in a dry place!

SERVICE CENTRE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified specialist or the service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan:

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (Waste Electrical and Electronic Equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead, it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

Disposal of batteries and power supplies

If the appliance contains an accumulator or a battery cell, these components contain substances which may be harmful to the environment when disintegrated. Before disposing of the appliance or replacing the accumulator or the battery cell, it is necessary to remove the original component and hand over to a suitable collection point.

CAUTION!

Accumulator or battery cell:

- **Do not throw into a fire!**
- **Never disassemble!**
- **Never shorten the wires!**

If leakage occurs when the accumulator is damaged, the electrolyte is of an alkaline character. If the electrolyte gets on your skin, wash it gently with soap and water and seek medical attention.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without prior notice and we reserve the right to change them.

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Warranty

The manufacturer (or importer) is responsible for ensuring the product complies with the requirements of applicable legal regulations as well as those of the relevant technical standards. Moreover, they are responsible for ensuring the product has the properties the manufacturer described in documents related to the goods or those reasonably expected by the customer with regard to the nature of the goods or based on advertising produced by the manufacturer, and further they are responsible for ensuring the product is fit for the purpose proposed by the manufacturer or that a product of the same type is normally used for.

The quality warranty term is 24 months from product takeover by the customer.

The warranty does not apply to wear and tear caused by regular use. The customer shall not be entitled to any warranty claims if, prior to taking the product over, they knew the product contained a defect or if the defect is attributable to the customer.

The warranty specifically does not apply:

- if the product installation, operation and service conditions stipulated in the product operating manual have not been adhered to,
- to malfunctions caused due to mechanical, heat or chemical damage, short circuit, over voltage or incorrect installation,
- to malfunctions caused by an inexperienced third-party intervention,
- to malfunctions caused by natural disaster,
- to malfunctions caused by insufficient or inappropriate maintenance in violation of the operating manual, including malfunctions caused by water and other sediments,
- to changes in colour of the heating surface or to scratching of the surface caused as a result of using the products in an unusual manner,
- to appearance and functional changes caused by exposure to sunlight, thermal radiation of water and other sediments,
- if the service life of certain product parts expires, e.g. for accumulators, bulbs, etc.

The warranty does not apply to any products and services provided along with the product (gifts, promotional articles, etc.).

Filing a complaint

A complaint against a product defect must be filed as soon as identified, yet no later than before the end of the warranty term.

The customer must file a product complaint with the dealer from which they have purchased the product, or with any authorised service centre, a list of which is included in the product package, or available at www.my-concept.com.

While filing a product complaint, the product must be duly cleaned and securely packed so as to prevent any damage during its transport to an authorised service centre, where relevant, unless the product is delivered in person.

The customer must submit proof of having concluded a purchase contract for the product by producing the receipt.

While filing their complaint, the customer must indicate the noted defect and identify the preferred complaint application method.

Complaints processing

As long as the noted defect may be removed, the user has the right to have the defect duly removed free of charge on a timely basis.

Where such a procedure is not reasonable with regard to the nature of the defect, the user may require to be supplied a new defect-free product (replacement), or, where the defect applies to a part of the product only, replacement of the part concerned. However, if replacement of the product or any part thereof is not proportionate with regard to the nature of the defect, especially if the defect can be removed without undue delay, the customer has the right to have the defect removed free of charge.

If the noted defect is not removable, or if the customer becomes entitled to replacement of the product or a part thereof, yet the replacement is not possible, for example due to the product having been sold out, the customer has the right to return the product (withdrawal from the contract).

The customer shall have the right to a new product (replacement) or to replacement of a part of the product even if the defect can be removed, provided they cannot properly use the product due to repeated occurrence of the defect or due to a high number of such defects. In such a case, the customer also has the right to return the product (by withdrawing from the contract).

If the product is not returned (the customer does not withdraw from the contract), or if the customer does not apply the right to a new defect-free product (replacement), to replacement of a part thereof or to repair of the product, they may request a reasonable discount. The customer also has a right to a reasonable discount if a new defect-free product cannot be supplied to them, or if a product part cannot be replaced or the product repaired unless the situation is remedied within a reasonable time limit, or if remedying the situation would create major discomfort on the part of the customer.

The seller, authorised service centre or a staff member authorised by them must decide about each complaint immediately or within three business days in complicated cases. This term does not include a reasonable period of time, depending on the type of product concerned,

required for the defect to be assessed by an expert. A complaint, including defect removal, must be dealt with without any undue delay, yet no later than within 30 calendar days of the complaint filing date, unless the seller and the customer agree on a later deadline.

On returning the product (withdrawing from the contract) the customer must return any accessories and documents supplied along with the product.

The customer does not have the right to keep the defective parts and components of the product replaced as part of a repair of the product.

This shall be without prejudice to any other rights the customer may have in relation to the purchase of the product.

Remark: Complaints against products damaged in transport are governed by the carrier's complaints procedure.

Manufacturer

Jindřich Valenta - Concept
 Vysokomýtská 1800, 565 01 Choceň
 Czech Republic
 tel.: +420 465 471 400
 fax +420 465 473 304
 Company ID No. 13216660
 email: servis@my-concept.cz
 www: www.my-concept.com

Product data

Model:	
Production number:	
Date of purchase:	Seal and signature of vendor:

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EU Declaration of conformity

Product / model Body composition smart scale VO4000, VO4001
Type: CS20N

Name and address of manufacturer: Jindřich Valenta – Concept
Vysokomýtská 1800
565 01 Choceň

This declaration of conformity is issued under the sole responsibility of the manufacturer.

Subject of the declaration: The product is intended for household. The product is designed for weighting the weight of people.

The product described above complies with the relevant harmonization legislation of the European Union:

directive 2014/35/EU (LVD),
directive 2014/30/EU (EMC),
directive 2011/65/EU (RoHS),
directive 2014/53/EU (RED)

Frequency range (GHz)	Max. Transmit power (dBm)
2,4 – 2,4835 GHz	≤ 20 dBm


Relevant harmonized standards used:

EN 62368-1:2014+A11:2017
EN 62479:2010
ETSI EN 301 489-1 V2.2.3
ETSI EN 301 489-17 V3.1.1
EN 55032:2015
EN 55035:2017
ETSI EN 300 328 V2.2.2

concept
Jindřich Valenta – Concept
Vysokomýtská 1800, 565 01 Choceň
Tel.: 465 471 400 - Fax: 465 473 304
IČ: 13216660 - DIČ: CZ6702171377

Signed and followed by the name:

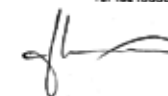In: Choceň 15.12.2020
Name and function: Dominik Jehlička, Product engineer SDA

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Verners VT Ltd.