HYUNDAI

OVET 739

NÁVOD K POUŽITÍ / NÁVOD NA POUŽITIE INSTRUKCJA OBSŁUGI / INSTRUCTION MANUAL HASZNÁLATÍ ÚTMUTATÓ

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
TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY. – Niebezpieczeństwo uduszenia. Tego worka nie należy

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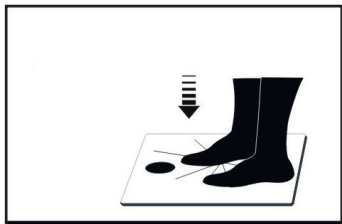
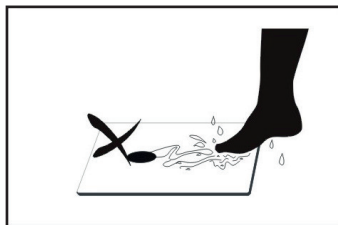
Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

I. SAFETY WARNING

- Before the first use, read the instructions for use carefully, look at the picture and store the instructions for use for the future. Consider the instructions for use as a part of the appliance and pass them on to any other user of the appliance.
- The appliance is intended for use in households only and for similar use! It is not intended for use in health-care facilities and for commercial use!
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never use the appliance if it malfunctions, if it fell to the ground and got damaged or if it fell into water. In these cases bring the appliance to a specialised service to have the safety and functions of the device checked.
- We recommend: Do not step onto the scale with wet feet or on wet surface of the scale, there is a risk of slipping.
- Never immerse the scale into water (even their parts)!
- Use the scale in a horizontal position only and at places where they cannot knock over, in sufficient distance from heat sources (e.g. stove, fireplace, heat radiators), sources of increased air humidity (e.g. sauna, bathroom, swimming pool) and equipment with strong electromagnetic field (e.g. microwave oven, radio, mobile phone).
- Unstable or soft surface under the scale can have negative effect on precision of weighing.
- In handling the scale, make sure injury of persons or damage of the appliance is prevented.
- Do not overload the appliance with weight higher than the maximum capacity! You will thus prevent damaging the appliance
- Do not use the fat analysis feature on this scale if you are using one of the following devices: pace maker, electrocardiograph and/or other medical device implanted in the body or device used as support.
- Do not put any objects on the surface of the scale during storing, they could be damaged.
- If you do not intend to use the scale for a long time, remove battery from them.
- If the battery is leaking, remove it immediately; it can damage the scale.
- Do not use worn and new batteries together and do not use rechargeable accumulators.
- Dispose of a discharged battery in a suitable way (see par. ENVIRONMENTAL PROTECTION).
- Never use the appliance for any other purpose than for that described in these instructions for use!

- Keep in mind that the data read from this scale are approximate and inaccurate in comparison to results of actual medicinal analyses. The exact determination of biometrical characteristics can be done only by a doctor via medical methods.
- Keep batteries and scale out of reach of children and incapacitated people. A person, who would swallow the batteries, must seek medical help immediately.
- Do not use during pregnancy.
- The producer is not responsible for damage caused by improper use of the appliance (e.g. breaking the board, etc.) and its guarantee does not cover the appliance in the case of failure to comply with the safety warnings above.

II. DESCRIPTION OF THE CONTROLS

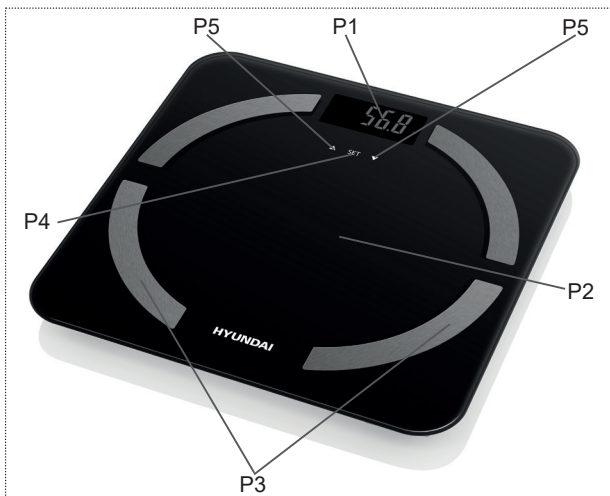- P1. Display
- P2. Glass surface
- P3. Electrodes
- P4. SET / ON button .
- P5. ▲ and ▼ buttons for setting
- P6. Batteries case cover

LCD SEGMENTUM / ICONS ILLUSTRATIONS:

Kg / lb / st:lb weight measurement units

P0 ...P9 number of memory

man

woman

MUS musculature contents (%)

BONE weight contents (%)

cm (ft: in) height measurement units

Age age

FAT body fat contents (%)

Underfat

Healthy

Overfat

Obese

TBW body water contents (%)

BMI body mass index

KCAL minimum energy consumption

SPECIFICATIONS:

Discrimination:	0.1 kg
Memory capacity:	10
Mass alignment:	from 60 to 250 cm
Age:	from 10 to 99 years
Minimum weight:	from 2.5 kg (5.5 lb / 0.4 st)
Maximum weight:	150 kg (330.6 lb / 23.6 st)
Body fat measuring unit:	0.1%
Body fat measuring range:	4 % – 50 %
Body water measuring range:	35 – 85 %
Muscle mass measuring:	10 % – 50 %
Bone measuring range:	2 % – 20 %

III. INSTRUCTIONS FOR USE**BEFORE USE**

- Open the battery cover on the bottom of the scale, insert the battery with the correct polarity and close the cover in the opposite direction.
- Place the scale on flat, stable and dry surface. Turn on the scale by tapping your foot on the weighing surface and check the weight unit displayed on the display (kg). If you wish to set unit to st or lb, adjust the buttons P5 on the scale kg ↔ st (st:lb) ↔ lb (or cm ↔ Ft;in).
- Place the scales on a firm, flat and even surface to ensure accuracy.

kg = kilogram

1 kg = 2.2 lb

st = stone (American weight unit)

1 st = 6.356 kg

lb = pound (British weight unit)

1 lb = 0.454 kg

Ft;in = feet, inches

1 FT = 30.48 cm; 1 In = 2.54 cm

WEIGHT MEASURING

You can use this scale only for weight measuring.

This is basic mode without body fat measuring.

- 1) Gently step on the scale, then the scale will be automatically turned on. Stand evenly on the scale without moving and wait until your weight shown on display is stable and locked. When you step off the scales it will automatically switch off after 6 seconds.
- 2) If no weighing is carried out while 0.0 is displayed the scales will automatically switch off.

BODY FAT MEASUREMENT

- 1) About Body Fat Percentage

Body fat percentage indicates body fat proportion in body, which is a criteria to judge one is overweight or not.

- 2) Testing Principle:

Impedance of fat is higher than any other body tissues such as muscle, when electrical signal of a certain frequency goes through body (BIA principle). Based on this principle, body fat percentage can be tested by sending an electrical signal through body.

BMI

The technology used in these scales, Bio-electrical Impedance Analysis (BIA) is more technologically advanced than the order Body Mass Index (BMI) method where fat % is calculated using only height and weight. BIA technology uses 5 parameters to get a more accurate resolution for the analysis (weight and biological resistance which are automatically measured by the scales together with the three personal parameters entered manually - sex, height and age). This index sets mutual ratio between the body weight and height. It can be used to determine whether the person suffers from underweight, overweight or obesity or whether his/her health condition is normal. Also with regards to heart and vascular diseases, body mass index is a relatively good risk index, particularly in healthy people. However, the index is for reference only. Experts recommend supplementing the calculation by at least measuring the waist, which makes determining disposition of the body fat more precise. BMI will not provide very reliable information for small children, elderly people and active sportsmen who have a lot of muscle mass as their results can be biased.

Minimum energy consumption CAL

It represents optimum amount of energy consumption necessary for maintaining the basic vital functions and balanced metabolism. The data on the display show the optimum daily consumption of energy that you give out at your weight, height, age and gender. You should also take in this amount of energy to keep your metabolism balanced. If you intend to reduce your weight, your daily energy intake should be a bit lower. However, remember to keep carbohydrates, fats and proteins balanced in your diet.

239 cal = 1 kJ

Important:

During body fat testing, please remove shoes and socks, etc., and clear your feet, in order to get the most accurate result. To ensure accuracy, please test at the same time of the day.

Test result may be inaccurate when the users are:

- after intense sport activity, overfeeding, or excessive dehydration,
- children under 10, or adults over 80,
- body building athletes, or other occupational athletes,
- pregnant women; Dialysis or edema patient,
- excessive drinking, eating, exercising, health conditions, used medicaments, women's menstruation cycle etc. influence the results of the analysis,
- figyeljen arra, hogy mérés közben ne legyen kontaktus a talpak, vádlik, vagy combok között. Ellenkező esetben lehetséges, hogy a mérés nem lesz pontos.

In order to ensure accurate results users should enter the required personal parameters; age, sex and height. These scale have the facility to store the personal parameters of 10 people.

SETTING AND SAVING PARAMETERS

Press **SET** to turn ON.

Press **▲** or **▼** to switch between Body Parameter P0 - P9. Press **SET** to confirm.

Press **▲** or **▼** to select male or female. Press **SET** to confirm.

Press **▲** or **▼** to indicate age. Press **SET** to confirm.

Press **▲** or **▼** to indicate height. Press **SET** to confirm.

Now all the personal data is saved (changed). The values are flashed and recorded.

MEASURING

- 1) Tap on the scales platform to switch on the product
- 2) Press **SET** P1 will flash. Press  or  to switch between Body Parameter P0 - P9.
- 3) Step onto the platform when the LCD display reads 0.0.
Important: For correct measurements and scale functions step on the scale barefoot.
- 4) First the weight is recorded, then „o o o o“ is displayed, indicating that the scales are analyzing. Please do not alight.
- 5) Once finished, the weight, body fat, body muscle, bone weight, BMI and KCAL will be displayed alternately for 2 times before the scales automatically turn off.
- 6) After this the scale is going turn off.

Compare the body fat reading against the <Body Fat Radio Chart> to confirm the body status.

Body Fat Ratio and Body water Ratio Chart

Age	Female		Male		Result
	Body fat (%)	Body water (%)	Body fat (%)	Body water (%)	
< = 30	4.0 - 20.5	66.0 - 54.7	4.0 - 15.5	66.0 - 58.1	Underfat
	20.6 - 25.0	54.6 - 51.6	15.6 - 20.0	58.0 - 55.0	Healthy
	25.1 - 30.5	51.5 - 47.8	20.1 - 24.5	54.9 - 51.9	Overfat
	30.6 - 50.0	47.7 - 34.4	24.6 - 50.0	51.8 - 34.4	Obese
< 30	4.0 - 25.0	66.0 - 51.6	4.0 - 19.5	66.0 - 55.4	Underfat
	25.1 - 30.0	51.5 - 48.1	19.6 - 24.0	55.3 - 52.3	Healthy
	30.1 - 35.0	48.0 - 44.7	24.1 - 28.5	52.2 - 49.2	Overfat
	35.1 - 50.0	44.6 - 34.4	28.6 - 50.0	49.1 - 34.4	Obese

Bone Weight Chart

Weight (kg)	Female			Male		
	< 45	45 - 60	> 60	< 60	60 - 75	> 75
Ref. bone weight (kg)	~ 3.0	~ 4.2	~ 6.5	~ 4.5	~ 6.0	~7.5

KCAL Chart

Age	Female		Male	
	Weight (kg)	KCAL	Weight (kg)	KCAL
6 - 17	50	1265	60	1620
18 - 29	55	1298	65	1560
30 - 49	60	1416	70	1680
50 - 69	60	1302	70	1561
70 - 100	50	1035	65	1397

BMI Chart

BMI	Underfat	Healthy	Overfat	Obese
	< 18.5	18.5 - 27	30 - 39	> 40

Due to the natural differences between people this chart should be used for reference purposes only. For the reason of differences between people, this table is only informative. Always consult your attending physician if any measuring results are abnormal. The results of measuring bones' weight and muscular proportions are simply tentative. Healthy man has approximately 40 % of muscular proportion and healthy woman has 30 % according to an average ratio. Women's weight of bones is among 12 % and men's weight of bones is among 15 % but it depends on many relevant factors.

ERROR AND OTHER LCD MESSAGE

LCD may show occasionally error message' Err ,. Simply turn off the unit for a brief moment, or remove and re-install battery as a total re-set. Then turn on again to resume operation. When the LCD screen displays:

- 1) „Lo“, you should install anew battery.
- 2) „Err 2“, it shows an error in the body fat percentage. Fat percentage exceeds the maximum capacity of 50 %, fat percentage is under 4 %.
- 3) „EEEE“ = The measured weight is higher than maximum allowed value of the scale.
- 4) „ErrH“ = Very high % of fat,
- 5) „ErrL“ = Very low % of fat.

If signs appear, such as if the display is flashing or if there is a faulty sign, try to move the scale away from the influence of a source of interference or ensure switching off the source for the time of using the scale.

CLEANING AND MAINTENANCE

- 1) Use a dry, soft fabric to clean the device. Or a lightly damp soft fabric with mild detergent to wipe clean the device. Do not use strong chemical solution such as alcohol, benzene for cleaning.
- 2) It could be slippery when the unit platform surface is wet, or semi-damp. Ensure the weighing platform surface is dry before using for safety reason.

TROUBLESHOOTING

The device is not working

- Check whether the scale is on a flat floor and they do not touch the wall.
- Check the rights position of battery and the polarity
- When warning indicator shows „Lo“ replace the battery.

IV. TECHNICAL SPECIFICATIONS

- Maximum weight capacity up to 150 kg
- Measurement accuracy 100 g
- Protection class III, CE Certification
- Size: 30 x 30 cm
- Weight (netto): 1.75 kg

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

Replacement of parts that require intervention in the electric part of the appliance must be carried out by a specialized service! Failure to comply with the instructions of the manufacturer will lead to expiration of the right to guarantee repair!

V. LEGISLATIVE & ECOLOGY

Disposal of old Electrical & Electronic Equipment (Applicable in the European Union and other European countries with separate collection systems).

This symbol on the product or on its packaging indicates that this product (or battery/accumulator) shall not be treated as household waste. Instead it shall be handed over to the applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The recycling of materials will help to conserve natural resources. For more detailed information about recycling of this product, please contact your local Civic Office, your household waste disposal service or the shop where you purchased the product.


WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY; DO NOT IMMERSE IN WATER OR OTHER LIQUIDS;


TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

The  symbol indicates a WARNING.

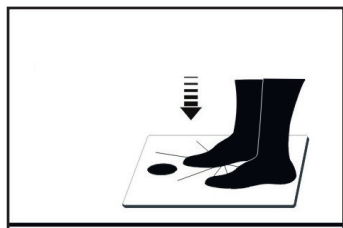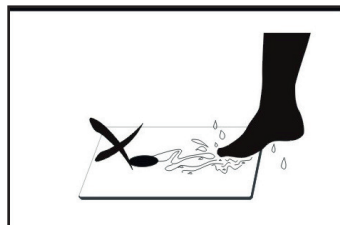AUTO STEP ON – just step on the scale and stand still, the scale will automatically turn ON and start taking measurements. Warning: Slippery when wet!

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Licensed by Hyundai Corporation, Korea.