SonicBrush GL UV

EN - USER MANUAL

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HR - UPUTE ZA UPOTREBU


EN - USER MANUAL

BOX CONTENTS

- Body of toothbrush
- 4x Standard brush heads
- Glass
- Charging station
- Travel case
- TrueLife SonicBrush UV Sterilizer
- Manual and safety instructions

PRODUCT DESCRIPTION


Sonic toothbrush with charging glass and travel case

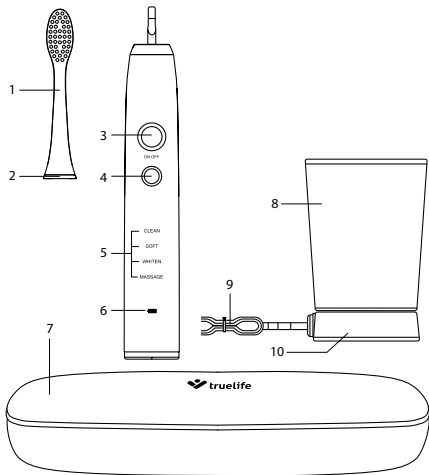- | | |
|--|-----------------------|
| 1. Replaceable brush head for sonic toothbrush | 6. Charging indicator |
| 2. Separating ring | 7. Travel case |
| 3. On/off button | 8. Glass |
| 4. Switch mode button | 9. Power cord |
| 5. Mode indicator | 10. Charging station |

UV Steriliser

1. Power cord input
2. Case cover
3. UV lamp
4. Button to unlock case cover
5. Holder for toothbrush/brush head

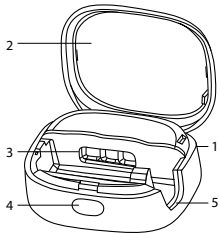INSTRUCTIONS FOR USE

Sonic toothbrush

1. Place the brush head onto the body of the toothbrush.
2. Select the cleaning mode using the button. If you intend to use the same mode as in the previous use, there is no need to reset the toothbrush.
 - **CLEAN** – Regular mode for brilliantly clean teeth
 - **SOFT** – Gentle brushing for sensitive teeth and gums (may also be suitable for users less accustomed to sonic cleaning)
 - **WHITEN** – For brushing with tooth polishing
 - **MASSAGE** – For gentle gum stimulation
3. Wet the brush head and apply toothpaste to it.
4. Place the bristles of the toothbrush on the teeth, near the gums.
5. Turn the toothbrush on using the on/off button.
6. Move the bristles back and forth over your teeth and slowly move between the individual teeth. At the same time, gently push the brush towards and between the teeth
7. For best results, divide the brushing into 4 sections. Take 30 seconds for each section. After a 30-second cycle, the brush will pause and alert you to move to the next section. After 2 minutes, the brush switches itself off.
8. Rinse the brush head under running water. Wipe the body of the brush with a damp cloth and then wipe it dry with a cloth or towel .
9. Put the body of the toothbrush back into the glass or store it in the travel case.

UV Steriliser

1. Use the button to open the lid of the UV steriliser .
2. Place the toothbrush or brush head in the UV steriliser. The bristles should face the UV lamp .
3. Close the lid. Sterilisation will begin automatically.
4. Sterilisation lasts 3 minutes. Do not open the lid during the sterilisation process.
5. Once the UV light switches off, you can open the lid and remove the toothbrush or brush head.

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