**Gebruiksaanwijzing
Mode d'emploi
Gebrauchsanleitung
Instructions for use
Instrucciones de uso
Manuale utente**

Alecto[®]

BC-09

NL	GEBRUIKSAANWIJZING	4
FR	MODE D'EMPLOI	10
D	BEDIENUNGSANLEITUNG.....	17
GB	USER MANUAL	24
ES	INSTRUCCIONES DE USO	30
IT	MANUALE UTENTE	36

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Introduction:

This ear thermometer is designed to measure your body temperature safely, quickly and accurately.

Body temperature can be measured more quickly and accurately in the ear than other places in the body. This means that the temperature measured with this ear thermometer is more accurate than the temperature measured in the armpit, mouth or rectum.

First use:

A battery has already been installed in the ear thermometer. Remove the insulating strip to activate the battery.

Power on / off:**Power on:**

Press button SCAN until all segments in the display are illuminated. Remove your finger from the button. The display will now show the last measured temperature.

Power off:

The thermometer will automatically switch itself off after one minute if no button is pressed again.

How to measure:

1. Switch the thermometer on as described.
2. Place the head of the thermometer in your ear (do not force and do not push too far in)
3. Press and hold button SCAN
4. You will hear a beep after a couple of seconds. You can now release the button and read the temperature on the display.

Display read-out:

 = temperature is being measured

 = temperature is below 34°C


 = temperature is below 38°C

 = temperature is above 38°C

 = temperature is above 43°C

 = last measured temperature

 = battery low

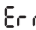 = ambient temperature is not between 16°C and 35°C:

Cleaning;


Temperature sensor:

Clean the temperature sensor with a cotton bud soaked in medicinal alcohol (96%, available from chemists) every time after use. NEVER clean the temperature sensor with water.

Thermometer:

You can clean the casing of the thermometer with a damp cloth. Do not use chemical cleaning agents.

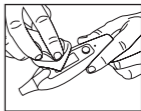Memory:

This ear thermometer remembers the last temperature measured. The thermometer will display this temperature when it is switched on.


Battery warning:

 : battery low, purchase a new battery (the thermometer can still be used in the meantime).

 + Lo : battery too low to function correctly (the battery should be replaced)

Battery replacement:

1. Remove the battery cover from the thermometer.
2. Use a pointed (non-metallic) object to remove the old battery and place the new battery with the + side up.
3. Place the battery cover back onto the thermometer.

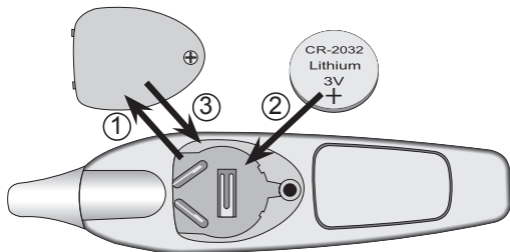Calibration:

The ear thermometer has been factory calibrated. Further calibration is unnecessary if the thermometer is used in accordance with the instructions. There is no means of readjusting the calibration on this thermometer.

Possible causes of the temperature measured being too low:

- The thermometer is not far enough or at an angle in the ear. Place the thermometer in the ear correctly and measure the temperature again.
- The temperature sensor is dirty. Clean it and measure the temperature again.
- The ear canal is blocked or dirty. Clean it or measure the temperature in the other ear.

For further questions and/or problems, contact the Alecto customer service.

Specifications:

Range: $34^{\circ}\text{C} \sim 43^{\circ}\text{C}$

Accuracy: $35.5^{\circ}\text{C} \sim 42^{\circ}\text{C}: \pm 0,2^{\circ}\text{C}$
 $<36^{\circ}\text{C} / >39^{\circ}\text{C}: \pm 0.3^{\circ}\text{C}$

Display resolution: $0,1^{\circ}\text{C}$

Ambient operating environment:

Temperature: $16^{\circ}\text{C} \sim 35^{\circ}\text{C}$

Relative humidity: up to 85%

Battery: CR2032, 3V, lithium

Battery lifetime: 2 year

3000 measurements

Environment:

Do not throw used batteries away. Hand them in at your local depot for domestic chemical waste.

At the end of this product's useful life, do not dispose of it with household rubbish, hand it in at a recycling collection point for electrical and electronic appliances.

Advice and warnings:

Never:

- Measure your temperature with the thermometer if you have an ear infection or other ear irritation.
- Measure your temperature within 30 minutes of physical activity or bathing.
- Disassemble the thermometer.
- Touch the thermometer sensor.
- Leave the thermometer in direct sunlight.

Always:

- Store the thermometer in a dry and dust-free environment and out of reach of children.
- Clean the thermometer well after use.

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**Wellkang Ltd