

Smart

Oral Irrigation
& Flossing Device

INSTRUCTION MANUAL



Welcome to



Please take a moment to read these instructions thoroughly so that you can enjoy all the features that make the Great Smile Store Smart Oral Irrigator our most powerful flosser yet.

The Great Smile Store Smart Oral Irrigator powers away plaque and debris between teeth and below the gum line to leave your mouth feeling incredibly fresh and clean.

Your Great Smile Store Smart Oral Irrigator is equipped with an advanced water pressure control system allowing you to personalize your irrigating and flossing experience. The Smart Oral Irrigator is rechargeable and has an easy-to-fill removable water tank.

Your Great Smile Starts Here.



Available Irrigating & Flossing Modes

The Great Smile Store Smart Oral Irrigator offers three different modes for varying oral care needs:

- **Normal** - Provides exceptional mouth cleaning and power to blast away plaque and debris
- **Soft** - Gentle, yet provides thorough cleaning for sensitive areas
- **Pulse** - Provides a gentle stimulation setting for your gums

Selecting Your Mode

- Power is activated by pressing the "Power" button (4).
- Mode selection is controlled by the "Mode" button (5). Pressing this button cycles through the working mode sequence Normal, Soft, Pulse.
- The Great Smile Store Smart Oral Irrigator automatically starts in the mode you last used.
- When you press the "Mode" button (5) the mode indicator will light up to show which mode you are using.
- To turn your Smart Oral Irrigator off before the 2 minute auto-timer finishes, press the "Power" button (4).




This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of this product in household waste at the end of its useful life. Disposal can take place at appropriate collection points provided in your area.



Two Minute Smart Auto Timer

The Great Smile Store Smart Oral Irrigator is programmed to run through a complete flossing & irrigation cycle in two minutes. The American Dental Association (ADA) recommends to brush twice a day and clean between your teeth with floss or oral irrigation once each day to remove plaque from all tooth surfaces. Plaque is responsible for both tooth decay and gum disease.

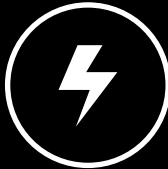

Connecting and Charging


It takes a minimum of 8 hours to fully charge the unit. A full charge last for approximately on month of regular use.

- Plug the adapter into the Smart Oral Irrigator..
- Plug the adapter into an electrical outlet.
- The "CH" indicator light on the unit will blink while the device is charging.



Using the Smart Oral Irrigator

- Attach nozzle to the unit.
- Open the water tank cap, hold the unit upright, fill the tank with water which temperature should not exceed 104°F, then close the cap. (Fig. 1)

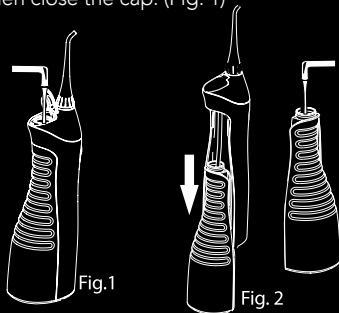


- Or, remove the water tank from the unit, hold the tank upright, fill the tank with water (temperature not is excess of 104°F), then attach the tank onto the unit. (Fig. 2)
- Press "Mode" key to select the mode you desire.
- Hold the unit in an upright position with nozzle tip pointed toward your teeth.
- Keep your mouth opened slightly so that the water can flow out.
- Press "Power" button to start irrigating.



After Using Smart Oral Irrigator

- Open the water tank cap and dispose of the remaining water in the tank.
- Press "Power" button to turn on the unit, and drain the remaining water from inside the appliance.
- Press "Power" button again to turn off the unit, and close the water tank cap.
- Wipe the unite dry with a soft cloth.

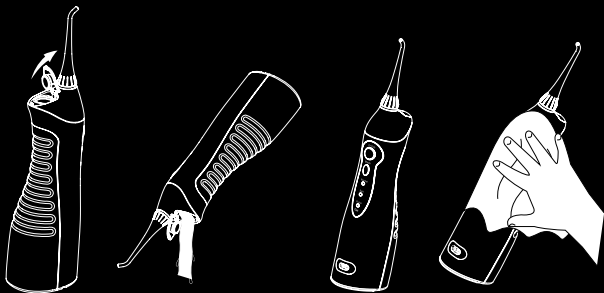

To prevent the growth of bacteria, make sure there is no remaining water inside the unit. If you're not going to use the Smart Oral Irrigator for some time, be sure wipe and dry the unit thoroughly before storing it.



Cleaning & Recommendations

Main Unit

- Remove the dirt by wiping the unit with a cloth and liquid soap.
- Do not submerge the main unit in water to clean.
- Wipe the unit dry with a soft cloth.

Water Tank Cap

- Wipe off any grime and dirt with a cloth.

Nozzle

- Wash with water and wipe with soft cloth.
- Replace the nozzle every 6 months.

Water suction tube

- Wash with water and wipe with soft cloth.
- Do not bend, pull or twist the water suction tube.

Tank

- Wash with water and wipe dry with soft cloth.
- Empty Dry the tank when not using for a long period.



IMPORTANT:

Use only water or neutral cleaning agent to clean the device. Do not wash using hot water of over 122°F. Keep the device away from high temperature area and avoid direct sunlight.

KEEP OUT OF REACH OF CHILDREN

Your Smart Oral Irrigator

- 1 Nozzle
- 2 Water Tank Cap
- 3 Smart Oral Irrigator Handle
- 4 Power Button
- 5 Mode Button
- 6 Mode Indicators
- 7 Water Tank
- 8 Charging Port

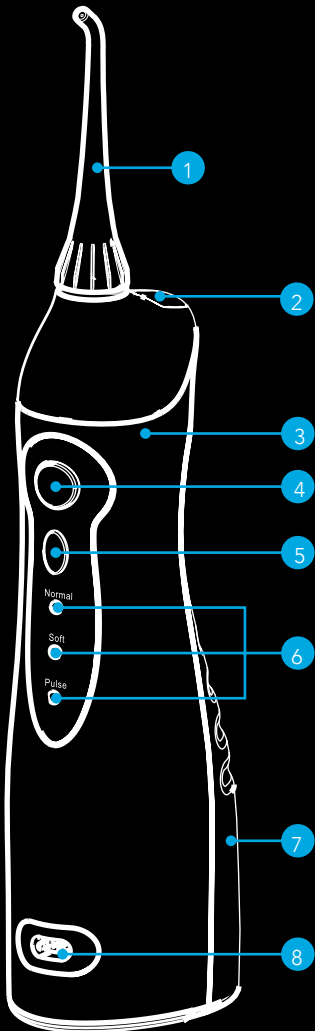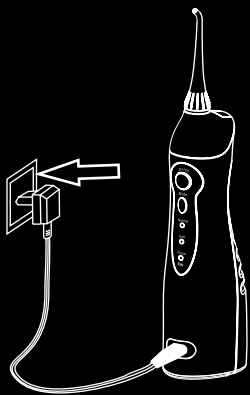



ONE - YEAR WARRANTY

This Great Smile Store product is warranted to be free of manufacturer's defects in materials or workmanship for one year from the date of purchase. Damage or wear resulting from an accident, misuse, abuse, commercial use, or unauthorized adjustment and/or repair are not covered under this warranty.

Should this product require warranty service (or replacement at our discretion) please contact client service to obtain a Return Merchandise Authorization number (RMA) and return instructions, proof of purchase is required. Products returned without a Great Smile Store generated RMA number will not be accepted and the sender will not receive a refund, replacement, or repaired product.

GREATSMILESTORE.COM

info@greatsmilestore.com

1.855.408.6669

There are no express warranties except as listed above. This warranty gives you specific legal rights and you may have other rights which vary from state to state.

DO NOT RETURN TO RETAILER, PLEASE CONTACT US DIRECTLY FOR WARRANTY SERVICE OR REPAIR, THANK YOU.

Thank You

FOR CHOOSING GREATSMILESTORE!

A healthy life should be full of wonderful surprises, so we would like to offer you

50% OFF

YOUR NEXT PURCHASE AT GREATSMILESTORE

YOUR SPECIAL DISCOUNT CODE:

[G S S M A R T 5 0]

Just enter your special discount code at checkout on your next purchase from:

 GREATSMILESTORE.COM

Copyright © 2015 Great Smile Store. All Rights Reserved.