 **wellneo**[®]

SONIC TOOTHBRUSH

HEALTHY TEETH AND GUMS

High frequency sonic vibrations

MODEL NUMBER: TB8061


User Manual

EN.....	3
AL	5
BG.....	7
BiH.....	11
CZ.....	14
EE	16
HR.....	18
HU	22
KS.....	25
KZ.....	27
LT	30


LV.....	32
MD.....	35
ME.....	37
MK	29
PL	42
RO.....	45
RS.....	47
RU	49
SI	52
SK	55
UA	57

(EN) PICTURES / (AL) FOTOGRAFIN / (BG) ИЛЮСТРАЦИИ / (BiH) SLIKE / (CZ) OBRÁZKY / (EE) JOONISED / (HR) SLIKE / (HU) KÉPEK / (KS) FOTOGRAFIN / (KZ) СУРЕТТЕР / (LT) PAVEIKSLIAUKAI / (LV) ATTĒLI / (MD) IMAGINI / (ME) SLIKE / (MK) СЛИКИ / (PL) RYSUNKI / (RO) IMAGINI / (RS) SLIKE / (RU) RYSUNKI / (SI) SLIKE / (SK) OBRÁZKY / (UA) МАЛЮНКИ


PICTURE 1

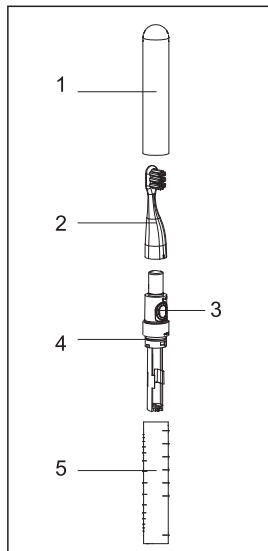PICTURE 2

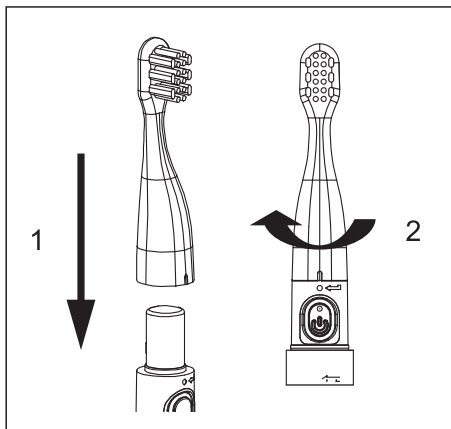PICTURE 3

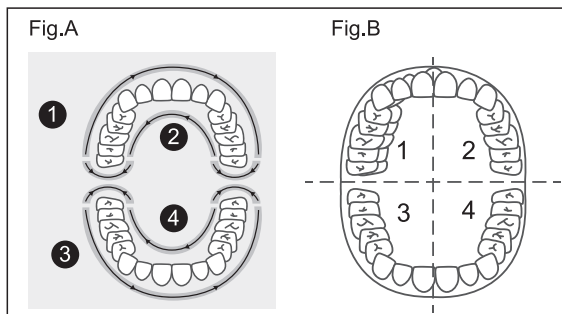Read this manual thoroughly before using and save it for future reference.

WARNINGS

1. This device is not a toy and should not be used by children or stored accessible to them.
2. This device is intended for household use only.
3. Only the batteries with a specification of 1.5V AAA can be used.
4. The instructions marked on the battery must be followed.
5. When replacing the battery, turn off the appliance first, then wipe out the water on the body. The battery must be installed correctly (polarity) and firmly.
6. Exhausted batteries are to be removed from the appliance and safely disposed of.
7. If the device is to be stored unused for a longer period (i.e. two weeks), the battery should be removed.
8. Do not use the battery which is out of date.
9. If the device is not working as it should, has received a sharp blow, has been dropped, damaged, left outdoors, or dropped into water, do not use it.
10. Do not attempt any modification or repair by yourself and ensure that any repair is conducted only by appropriately qualified technician.
11. In case of malfunction and/or damage do not attempt any repair by yourself and ensure that the product is properly disposed in accordance with local laws.
12. Do not immerse the toothbrush into the high-temperature water.
13. This device can be used by children aged from 8 years and above as well as by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved. Children should not play with the device. Cleaning and user maintenance should not be done by children unless they are older than 8 and supervised.
14. Do not press the tooth or gum by the toothbrush with excessive force and avoid brushing the same position for long time to avoid any injury.
15. Please consult the doctor whether you can use the toothbrush.
16. Oral bleeding may occur when using the toothbrush for the first time or improperly. If it lasts for more than a week, please consult the doctor for suggestion.

17. Please consult the doctor for suggestion before using the toothbrush if your body is equipped with heart pacemaker or other electronic devices.
18. Stop using a brush head with missing or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
19. Although the device has been checked its usage and consequences are strictly user responsibility.

PARTS (See Picture 1)

1. Protective bristle cover
2. Brush head
3. Power button
4. Main body
5. Bottom cover

BEFORE FIRST USE

Installing the battery

1. Detach the bottom cover by turning it counter clockwise and down.
2. Install the battery by turning it according to correct polarity and inserting it in and down.
3. Reposition the bottom cover by aligning the marks first and then rotating it clockwise.

Installing/replacing the brush head (See Picture 2)

1. Align the marks on brush and main body before sliding the brush head onto the main body.
2. Turn the brush head clockwise until it aligns with the mark on main body.

OPERATION

Wellneo® Sonic Toothbrush is a premium quality brush which uses DuPont registered Tynex® nylon bristles to suit teeth topography for thorough cleaning of teeth and gums and hard to reach areas.

This toothbrush has 2 different performing modes:

1. CLEAN- daily/ regular teeth brushing and removing of plaque
2. SOFT-sensitive cleaning and massage (gum care)

NOTE: Beginners may experience a tickling/tingling sensation but normally get used to it within a week.

BRUSHING INSTRUCTIONS

Follow these simple steps to maximize your brushing experience.

1. Remove the protective bristle cover.
2. Moisten the brush under running water.
3. Put some toothpaste to the moistened brush and place it in your mouth.

NOTE: Wait until the brush is in the mouth before switching it on, to avoid spraying.

4. Place the toothbrush bristles against the teeth at a slight angle towards the gumline.
5. Press the power on/off button to switch ON the Wellneo® Sonic Toothbrush. It will start operating in "CLEAN" mode with high vibrations. You can switch to "SOFT" mode with lower vibrations if power button is pressed again. Choose the mode which suit you.
6. Apply light pressure to maximize Wellneo® Sonic Toothbrush's effectiveness and let the toothbrush do the brushing for you. **Do not scrub.**
7. Gently move the brush head slowly across the teeth in a small back and forth motion so the bristles reach between your teeth. Continue this motion throughout your brushing cycle (see Picture 3).

NOTE: To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections (see Picture 3).

8. For optimum results we suggest one of the below cleaning methods:

- a. See Picture 3, Figure A: Begin brushing in section 1 (outside top teeth) and brush for 30 seconds before moving to section 2 (inside top teeth). Continue brushing in section 3 (outside bottom teeth) and brush for 30 seconds before moving to section 4 (inside bottom teeth). All 4 sections are brushed in 2 minutes.
- b. See Picture 3, Figure B: Begin brushing in section 1 (upper right quadrant) and brush for 30 seconds before moving to section 2 (upper left quadrant). Continue brushing in section 3 (bottom right quadrant) and brush for 30 seconds before moving to section 4 (bottom left quadrant). All 4 sections are brushed in 2 minutes.
9. Toothbrush turns OFF automatically after 2 minutes. Rinse the brush head with running clean water and shake it well to remove the water, dry the main body and bottom cover with a towel, finally, close with the protective bristle cover.

AUTO-TIMER

Dentists recommend brushing for 2 minutes at least twice a day, 30 seconds per quadrant. Our timer notifies you each 30 seconds when it is time to move to the next quadrant and shuts off after two minutes so you give your teeth the optimal brushing! Brushing cycle is complete when it automatically turns OFF the toothbrush. To clean your mouth more completely, you may also brush your tongue, with the toothbrush turned ON or OFF, as you prefer.

CLEANING AND MAINTENANCE

Before cleaning, switch off the device and make sure that the bottom cover has been closed well.

Rinse the upper cover, handle and brush head with clean water (with temperature below 45°C).

CAUTION: Do not immerse the toothbrush into water. Never use a poisonous cleaner. Do not detach the bottom cover for or while cleaning.

CAUTION: Do not use chemical, steel, wooden or abrasive cleaners to clean the toothbrush.

TIP: Replace the brush head every 3 months or sooner if signs of wear appear.

TROUBLESHOOTING GUIDE:

Problem	Reason	Solution
Product doesn't work.	The batteries have ran out.	Change the batteries.
	The battery installation is not correct.	Please check the position of the batteries.

If any type of other problem occurs, please contact your local supplier.

TECHNICAL DATA

Product must be kept dry and not exposed to direct sunlight during transport and storage.

Operating temperature and humidity: 0°C-50°C, 30%-90% RH

Storage temperature and humidity: -10°C-60°C, 30%-90%

Service life: 1 year

Product name: WELLNEO SONIC TOOTHBRUSH / ID: 110002937 / Model no.: TB8061

Voltage: 1.5 V-1x AAA battery (not included) / Power: 0.3 W / Current: 0.15 A / Noise: <65 dB

Country of origin: P.R.C.

Trademark owner/Importer to EU/ Distributor for Studio Moderna: Top Shop International S.A., Via Ferruccio Pelli 13, 6900 Lugano, Switzerland

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Production date and serial number can be found on the product.

Wellneo is a registered trademark of Top Shop International S.A.

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