exacta®

DELUXE

BEDIENUNGSANLEITUNG
OPERATING INSTRUCTIONS
MODE D'EMPLOI
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
INSTRUCCIONES DE MANEJO
MANUAL DE INSTRUÇÕES
BRUKSANVISNING
BRUGSANVISNING
KÄYTTÖOHJEET
KEZELÉSI ÚTMUTATÓ
INSTRUKCJA OBSŁUGI
NÁVOD K POUŽITÍ
ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ
KULLANMA KILAVUZU
ΟΔΗΓΙΕΣ ΛΕΙΤΟΥΡΓΙΑΣ
NAVODILA ZA UPORABO
UPUTA ZA UPORABU




Quality by **SOEHNLE**

- DE** Nicht geeignet für Personen mit elektronischen Implantaten (z.B. Herzschrittmacher).
- EN** It is not suitable for persons with electronic implants (pacemaker, etc.).
- FR** Les personnes qui portent des implants médicaux (p. ex. un stimulateur cardiaque etc.) ne doivent pas utiliser cet appareil.
- IT** Non è adatto per persone con impianti elettronici (pacemaker, ecc.).
- NL** Niet geschikt voor personen met elektronische implantaten (pacemaker, enz.).
- ES** No es apropiado para personas con implantes electrónicos (marcapasos, etc.).
- PT** Não é adequado para pessoas com implantes electrónicos (pacemakers, etc.).
- DA** Ikke egnet til folk med elektroniske implantater (pacemaker, etc.).
- SV** Ej lämplig för personer med elektroniska implantat (Pacemaker etc).
- FIN** Tuote ei soveltu henkilöille, joilla on elektronisia istutteita (sydämen tahdistin tms.).
- HU** Nem alkalmazható olyan személyeknél, akielektronikus implantátumokkal (például szívritmusszabályozóval stb.) rendelkeznek.
- PL** Nie nadaje się dla osób, którym wszczepiono elektroniczne implanty (np. rozruszniki serca).
- CS** Není vhodný pro osoby s elektronickými implantáty (kardiostimulátory, atd.).
- RUS** Никогда не используйте данный прибор вместе с электронными медицинскими приборами, например, кардиостимуляторами.
- TR** Elektronik parçalar (kalp pili, vs.) taşıyan insanlar için uygun değildir.
- GR** Δεν είναι κατάλληλο για άτομα με ηλεκτρονικά εμφυτεύματα (βηματοδότης, κ.λπ.).
- SLO** Ni primerno za osebe z elektronskimi implantacijami (spodbujevalec srca itd.)
- HR** Nije pogodan za lica sa elektroničnim implantatima (bajpas, itd.).

- DE** Achtung! Rutschgefahr bei nasser Oberfläche.
- EN** Attention! Slipping hazard in case of wet surface.
- FR** Attention ! Risque de glissade en cas de surface humide.
- IT** Attenzione! Pericolo di scivolare in caso di superficie bagnata.
- NL** Opgelet! Slipgevaar op nat oppervlak.
- ES** ¡Cuidado! Peligro de resbalar encima de una superficie mojada.
- PT** Atenção! Perigo de escorregar em superfícies molhadas.
- DA** Pas på! Fare for at glide på våd overflade.
- SV** Observera! Risk att halka på våta ytor.
- FIN** Huom! Liukastumisvaara märällä pinnalla.
- HU** Figyelem! Ha a felület nedves, fennáll az elcsúszás veszélye.
- PL** Pozor! V případě vlhkého povrchu nebezpečí uklouznutí.
- CS** Uwaga! Niebezpieczeństwo poślizgnięcia się na mokrej powierzchni.
- RUS** Внимание! Опасность одскальзывания на мокрой поверхности.
- TR** Dikkat! Islak zeminlerde kayma tehlikesi vardır.
- GR** Προσοχή! Κίνδυνος ολίσθησης σε υγρές επιφάνειες.
- SLO** Pozor! Nevarnost zdrsa na vlažni površini.
- HR** Pozor! Opasnost od klizanja ako je površina mokra.

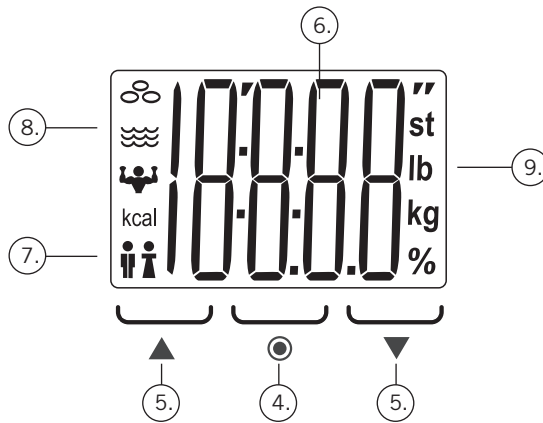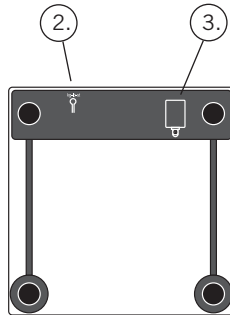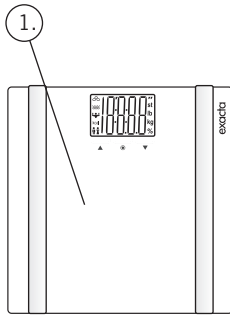DE

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EN

Entry of personal data

1. Push \odot key.
2. Use $\blacktriangle/\blacktriangledown$ keys to select memory slot and confirm (\odot).
3. Then use $\blacktriangle/\blacktriangledown$ keys enter body height. Use \odot to confirm.
4. In the next step enter the age and confirm.
5. Finally select sex and confirm.
6. Place scale on hard, level surface and wait until the scale switches off automatically.

Note: Changing or deleting data is done by overwriting existing data.

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EN

Guarantee

For 3 years as of the date of purchase, Exacta guarantees the removal of defects due to material or processing errors, free of charge, by way of repair or replacement. In case of a guarantee claim, please return your scale to your seller, together with purchase receipt and guarantee card.

EC Conformity

This device complies with the applicable EC Directive 2004/108/EC.

Battery disposal EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries containing hazardous materials are marked with these symbols:

- Pb = Battery contains lead
- Cd = Battery contains cadmium
- Hg = Battery contains mercury

Disposal of electric and electronic devices

EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

Consumer service

If you have questions or comments, the following points of contact are available:

Germany	Ph: (08 00) 5 34 34 34
International	Ph.: +49 26 04 97 70
Monday to Friday	08:30 to 12:00 am

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Leifheit AG